Bathroom Scale

Art.no 44-1845, -1, -2 Model BG198

Please read the entire instruction manual before use and save it for future reference. We reserve the right for any errors in text or images and for making any necessary technical changes to this document. If you should have any questions concerning technical problems please contact our Customer Services.

Safety

- This product is a precision instrument and should be handled with care and protected from blows and rough treatment.
- Place the scale in the same place every time you weigh yourself.
- For accurate readings, the scale should be placed on a hard, level surface. An uneven floor beneath the scale can cause the readings to fluctuate.
- Weigh yourself without clothing at the same time of day, preferably in the morning before breakfast. Place your feet in the same place every time you weigh yourself.
- The displayed readings are rounded up/down to the nearest 100 g.
- The scale can be cleaned using a damp cloth, if and when required. Make sure that no moisture gets inside the product. Never use solvents or strong detergents.
- Keep the scale away from water and splashes; otherwise the electronic components can be damaged.
- Do not store the scale standing on its edge as the display will be on all the time and drain the battery.
- Remove the battery from the scale if it is not to be used for an extended period.
- The scale can be affected by changes in temperature or by being moved. Inaccurate readings may occur if the scale is too cold when used. Allow the scale to stand and acclimatise for approx. 24 hours before using it for the first time.


Product description

Digital glass topped bathroom scale.

Battery

- Remove the isolating plastic tab between the battery and the battery holder before first use.
- Changing the battery: Carefully press the locking tab in to release the battery cover and then remove the battery. Insert a new battery (CR2032). Refit the battery cover.
- When the battery needs changing, “Lo” will appear on the display.

Buttons and functions

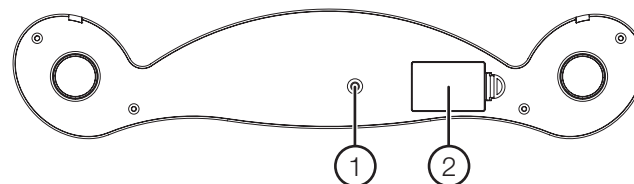1. [UNIT] Unit of measurement selector
2. Battery cover

Operation

1. Press [UNIT] to select the unit of measurement.
2. The scale should be placed on a level surface.
3. Stand carefully on the scale and stand still to allow your weight to stabilize.
4. Your measured weight will show in the display.
5. The scale will turn off automatically after 10 seconds.

If the scale is overloaded “Err” will appear on the display.

Care and maintenance

Clean the product by wiping it with a soft, damp cloth. Never clean the product using solvents or abrasive cleaning agents.

Responsible disposal

This symbol indicates that this product should not be disposed of with general household waste. This applies throughout the entire EU. In order to prevent any harm to the environment or health hazards caused by incorrect waste disposal, the product must be handed in for recycling so that the material can be disposed of in a responsible manner. When recycling your product, take it to your local collection facility or contact the place of purchase. They will ensure that the product is disposed of in an environmentally sound manner.

Specifications

Weight capacity	150 kg/330 lb
Graduation	100 g/0.2 lb
Batteries	1 × CR2032 (included)

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