Mini-Massager

Art.no 44-2720 Model M20300AD00-EU

Please read the entire instruction manual before using the product and then save it for future reference. We reserve the right for any errors in text or images and any necessary changes made to technical data. In the event of technical problems or other queries, please contact our Customer Services.

Safety

Note: The product is NOT meant to be used:

- On animals or children.
- On burn injuries, open cuts or sores.
- On pregnant women, recently operated or sensitive body parts such as the head and joints.
- On people bearing pacemakers, implants or other physical aids, or if you are unsure whether the product is appropriate for you. In these cases, please consult a doctor before using.
- The product is only intended to be used on individual body parts on the human body.
- This product must not be used by anyone (including children) suffering from physical or mental impairment. This product should not be used by anyone that has not read the instruction manual unless they have been instructed in its use by someone who will take responsibility for their safety.
- Never let children play with the product.
- The product is only intended for personal domestic use and only in the manner described in this instruction manual.
- Never use accessories not approved by the manufacturer.
- Never attempt to modify the product in any way.
- Never use the product if it is damaged or functions abnormally.
- Keep the product away from hot surfaces.
- The product must not be immersed in water or other liquid.

Health considerations


- Contact a doctor before using this product if you are pregnant, use a pacemaker or have any of the following medical diagnoses: Diabetes, heart disease, malignant melanoma or a condition which requires substantial rest.
- If soreness, rashes or pain result after use, discontinue use immediately and contact your doctor.
- Never use the product on swollen or inflamed areas. Do not use on areas where there are changes in the skin or on open sores.
- Do not use the massager more than 15 min during any session.

Battery safety

- Always use the correct type batteries, 2 × AAA/LR03. Ensure that the polarity is correct before inserting them.
- Do not mix different types of batteries, e.g. rechargeable batteries w/non-rechargeable or old with new.
- Remove the batteries if the appliance is not to be used for a long time.
- Do not throw batteries onto the fire, they can explode.

Operating instructions

1. Turn the bottom section so that its arrow is aligned with the arrow on the housing.
2. Remove the bottom section.
Note: The waterproof seal which allows the massager to be used in the shower may make removing the bottom section somewhat difficult.
3. Remove the battery cover and insert 2 × AAA/LR03 batteries. Follow the markings on the battery cover to ensure correct polarity.
4. Refit the battery cover.
5. Refit the bottom section and twist it to lock it securely in place.

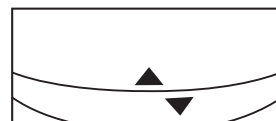Care and maintenance

- Clean the mini-massager by wiping it with a damp cloth. Use only mild cleaning agents, never solvents or corrosive chemicals.
- Whenever the mini-massager is not to be used for an extended period, the batteries should be removed and stored in a dry, dust-free environment, out of children's reach.

Responsible disposal

This symbol indicates that this product should not be disposed of with general household waste. This applies throughout the entire EU. In order to prevent any harm to the environment or health hazards caused by incorrect waste disposal, the product must be handed in for recycling so that the material can be disposed of in a responsible manner. When recycling your product, take it to your local collection facility or contact the place of purchase. They will ensure that the product is disposed of in an environmentally sound manner.

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