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Art.no 34-1019 Model M101B


Please read the entire instruction manual before using the product and save it for future reference. We reserve the right for any errors in text or images and any necessary changes made to technical data. If you have any questions regarding technical problems please contact Customer Services.

## Safety

- Never use the product for anything other than its intended use.
- Never place or store the product where there is a risk that it can fall into liquid e.g. the bath or washbasin.
- Never use the product as a toy.
- Always turn the product off when it is not in use.
- Never use the product if it is damaged or does not function correctly.
- Keep the product away from sources of heat.
- This product is intended for private use only.
- This product must not be used by anyone (including children) suffering from physical or mental impairment. This product should not be used by anyone without sufficient experience or knowledge unless they have been instructed in its use by someone who will take full responsibility for their safety. Never let children play with the product.
- This product is not intended for commercial use. Do not use this product as a substitute for medical expertise.

## Battery

- Use only the type of battery specified.
- Insert the batteries according to the polarity markings in the battery compartment. Refer to the markings in the battery compartment.
- Do not mix different types and qualities of batteries e.g. old with new.
- If the product is not to be used for an extended period, remove the batteries to avoid any possible risk of leakage.
- When the batteries need changing, open the battery compartment lid by turning it anticlockwise.
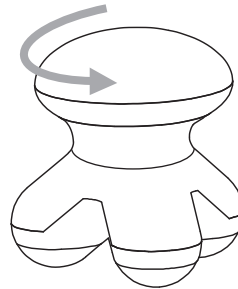
## Use

- Press and hold the handle down to turn the massager on.
- Do not use the massager for more than 15 minutes at a time.
- Turn the massager off by releasing the handle.

## NB

If you have a pacemaker, are pregnant or have diabetes, you should consult a doctor before using the massager.

## Tip

Drink a few glasses of water before and after a massage since this assists in extracting waste from the body.

## Care and maintenance

Clean the massager every now and then with a slightly moistened cloth. Use a mild detergent and avoid solvents and strong, abrasive cleaning agents.

## Disposal

Follow local ordinances when disposing of this product. If you are unsure of how to dispose of this product, please contact your local authority.

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