CAPERRE

44-9770

RED LIGHT
RECOVERY BELT

Red light recovery belt

Art.no 44-9770

Model HML-BL-DC-SW

Please read the entire instruction manual before using the product and then keep it for future reference. We accept no liability for any errors in the text or images and we reserve the right to make any necessary changes to technical data. If you have technical problems or other queries, please contact our customer service team.

Safety

- The product must not be used by people who are not sensitive to heat or are vulnerable in other ways, as they cannot react to overheating (for example, diabetics, people with skin changes caused by disease, people with scars in the treatment area or people who have taken medication or drunk alcohol).
- The product can be used by children aged 8 and over and by people with any form of disability or with a lack of experience and knowledge which could put their safety at risk, provided that they have been given instruction in using the product safely and understand the potential risks. Never let children play with the product. Children must not be allowed to clean or maintain the product unless they are under adult supervision.
- Never let children play with the product.
- The product can be cleaned or maintained by children aged 8 or over if they are under adult supervision.
- This product is intended for private domestic use.
- Do not use the product when it is folded or rolled up.

ENGLISH

- **Warning:** Do not use the product near baths, showers, washbasins or other items that contain water or other liquids.
- Consult a doctor if you have a pacemaker, skin disease or similar symptoms.
- Do not pull, twist, bend or fold the cables.
- Never use the product if it is damaged or not functioning correctly.
- If the mains lead is damaged, it must be replaced by a lead with the same specifications as the original.
- Do not use the product with other heat sources, such as hot water bottles, heated cushions or similar items.
- Examine the product regularly to detect any signs of wear or damage. If the product is damaged in any way, it must not be used.
- The product is not intended for medical use in hospitals.
- The product is not suitable for use with direct skin contact by anyone who matches any of the following descriptions:
 - Heart disease
 - Mental illness
 - Abnormal blood pressure
 - Sunburnt skin
 - Keloid scars
 - Oversensitivity to light
 - Discolouration and contact dermatitis
 - Infectious diseases
 - Immunological diseases
 - Diabetes
 - Melanomas on the face
 - Abnormal wound healing
 - Malignant tumours

ENGLISH

- People who are taking medication, cancer patients, especially with skin cancer or precancerous moles
- Pregnant women, breastfeeding women and other sensitive people
- Patients with capillaries caused by long-term use of steroid hormone medications or by liver dysfunction

Responsible disposal

This symbol indicates that the product should not be disposed of with general household waste.

This applies throughout the entire EU. To prevent any harm to the environment or health hazards caused by incorrect waste disposal, the product must be handed in for recycling so that the materials can be disposed of in a responsible way.

To recycle the product, take it to your local recycling facility or contact the retailer. They will ensure that the product is disposed of in a way that does not harm the environment.


Specifications

Input	5 V DC/2 A, max 10 W
Connection	USB-A
Wavelength	Red: 660 nm Near infra red: 850 nm
Dimensions	1290 × 220 × 8 mm
Weight	525 g

Product description

1. Hook and eye closure
2. Light panel
3. [⏻] On/Off
4. USB-A connector
5. Pocket for portable power source

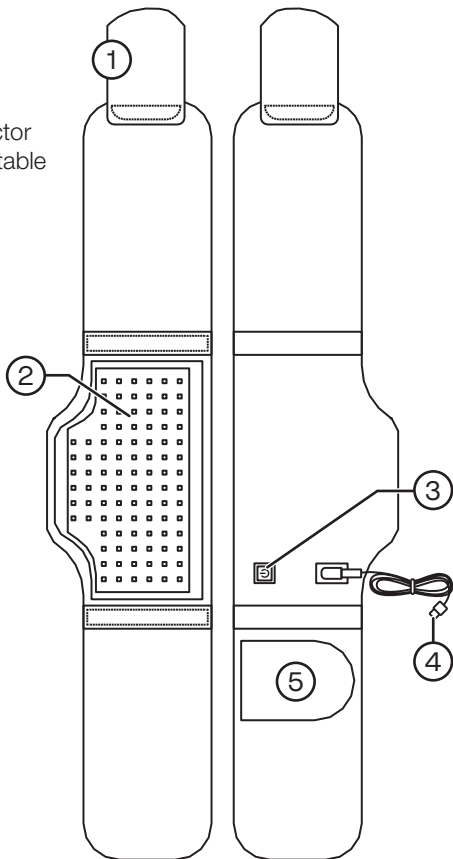Use

1. Fasten the belt around your body with the light panel either against your stomach or against your back, depending on which area you want to treat.
2. Connect the belt's USB-A connector (4) to the desired power source.
3. If a portable power source is used, such as a power bank, it can be conveniently placed in the pocket (5).
4. Hold [⏻] (3) for 3 seconds to turn on. It starts at the highest intensity.
5. Press [⏻] repeatedly to select the intensity as follows:
High → Medium → Low ↻
6. The belt has a timer that turns it off after 30 minutes. To turn it off manually, press and hold [⏻] for 3 seconds.

Cleaning and maintenance

- Wipe the light panel with a soft, slightly damp cloth. Never use solvents, abrasive cleaners or corrosive chemicals.
- Never immerse the product in water or rinse it under a tap. It is not waterproof.

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Post Clas Ohlson AB, 793 85 INSJÖN
Osoite Clas Ohlson Oy, Kaivokatu 10 A,