CAPERE

44-9569

LIGHT THERAPY
MASK

Soft light therapy mask

Art. no. 44-9569

Model G82


Please read the entire instruction manual before use and save it for future reference. We reserve the right for any errors in text or images and for making any necessary technical changes to this document. If you should have any questions concerning technical problems please contact our Customer Services (see address details on the back).

Safety

- This product is intended for use on the face only.
- The product is not intended for use as a medical treatment.
- The product must not be used by children under the age of 12.
- The product can be used by children aged between 13 and 18 under the supervision of an adult.
- Never use the product in the bathroom or in other damp environments. Never immerse the product in water or other liquid.
- Never use the product if it is not functioning correctly or is damaged in any way. Do not use the product if you have dropped it on the floor or if it has fallen into water or other liquid.
- Never stand or sit on the product.
- Never put the product close to a heat source, such as a cooker, woodburning stove, oven or open fire.
- Never cover the product while it is in use.
- Never try to open the housing or repair or modify the product in any way.
- **Risk for eye damage.** Do not look directly at the LEDs.

- The product is not suitable for use by anyone who has any of the following conditions or who matches any of the following descriptions:
 - Silicone allergy
 - Heart disease
 - Mental illness
 - Abnormal blood pressure
 - Sunburnt skin
 - Keloid scars
 - Oversensitivity to light
 - Discolouration and contact dermatitis
 - Infectious diseases
 - Immunological diseases
 - Diabetes
 - Melanomas on the face
 - Abnormal wound healing
 - Malignant tumours
 - People who are taking medication, cancer patients, especially with skin cancer or precancerous moles
 - People who have had botox injections in the face within the last three months
 - Pregnant women, breastfeeding women and other sensitive people
 - Patients with capillaries caused by long-term use of steroid hormone medications or by liver dysfunction
- If you have permanent or temporary fillers in the area that will be covered, consult a doctor before using the product.
- If you are undergoing hospital treatment, consult a doctor before using the product.
- A certain amount of redness can occur during the treatment. This is normal and should disappear within 1–2 hours. If you feel a burning sensation or other discomfort during the treatment, stop the treatment immediately.

Product description

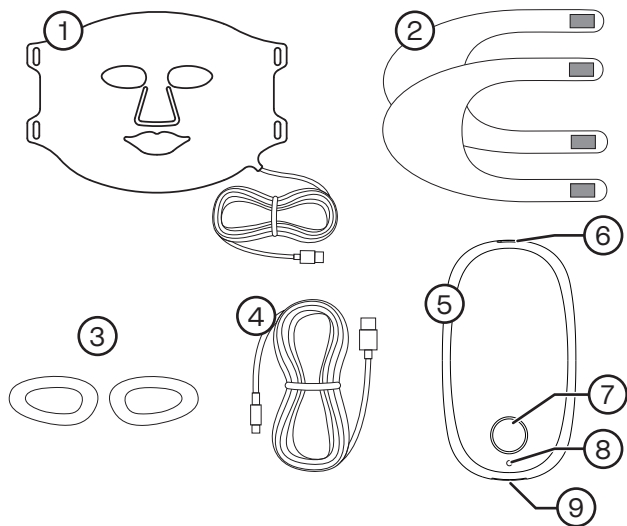- | | |
|----------------------------------|--------------------------|
| 1. Light therapy mask | 5. Control box |
| 2. Strap | 6. USB-C input port |
| 3. Eye cups for light protection | 7. Multi-function button |
| 4. USB-A to USB-C cable | 8. Indicator |
| | 9. USB-A output port |

Charging

Charge the product fully before using it for the first time.

1. Connect a USB cable to the USB-C port (6) on the control box and to a suitable USB charger.
2. The LED indicator (8) will pulse green while the battery is charging. When the battery is fully charged, the LED will show a steady green light.
3. The LED will flash red when the battery needs charging.

Photosensitivity test:

Note You should carry out a photosensitivity test before using the product for the first time:

Switch the product on and put it on the inside of your forearm. Leave it there for 10 minutes. If any type of redness or rash appears, we advise against using the product, as your skin is probably extra sensitive. Consult a skin specialist before continuing to use the product.

Using the product

Note The product must only be used with the control box provided. The product will be unsafe if you use it in any other way.

Before using the product for the first time: Attach the strap to the light therapy mask, fit the eye shells in the mask's eye holes and connect the mask's USB-A port to the USB-A output port on the control box.

1. Cleanse your face. Skincare products or makeup can contain ingredients that have a negative reaction to light.

2. Put the mask on your face and fasten the strap with Velcro fastener. Adjust the mask so that it fits firmly but comfortably. **Note:** Make sure that the eye cups fit snugly around your eyes. If they do not, keep your eyes closed during treatment.
3. Press and hold the multi-function button (7) for 3 seconds to start the mask. The mask starts with multicoloured light.
4. Press the multi-function button repeatedly to toggle between the following light settings:
Multicoloured light (850 nm) → Red light (660 nm) → Blue light (415 nm) → Yellow light (585 nm) ↺
5. The mask has a timer function that switches it off automatically after 10 minutes. To switch off the mask manually, press and hold the multi-function button for 3 seconds.
6. You can follow your usual skincare routine after a treatment. **Note** A certain amount of redness can occur during treatment. This is normal and should disappear within 1–2 hours.

Cleaning and maintenance

- Wipe the mask with a soft, slightly damp cloth. Never use solvents, abrasive cleaners or corrosive chemicals. Be careful not to damage the light panel or silicone surface.
- Never immerse the mask in water or rinse it under a tap. It is not waterproof.

Responsible disposal

This symbol indicates that the product should not be disposed of with general household waste.

This applies throughout the entire EU. To prevent any harm to the environment or health hazards caused by incorrect waste disposal, the product must be handed in for recycling so that the materials can be disposed of in a responsible way.

To recycle the product, take it to your local recycling facility or contact the retailer. They will ensure that the product is disposed of in a way that does not harm the environment.

Specifications

Input	5 V DC/2 A
Charging port	USB-C
Charging cable	USB-A – USB-C
Battery	3000 mAh Li-po
Charging time	approx 2 h
Battery life	4–5 uses of 10 minutes each
Wavelength	Multi: 850 nm Red: 660 nm Blue: 415 nm Yellow: 585 nm
Dimensions	311 × 217 × 5 mm
Weight	225 g

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