Massager with heat treatment

Art. no. 44-9672 Model FE-7111

Please read the entire instruction manual before using the product and then keep it for future reference. We accept no liability for any errors in the text or images and we reserve the right to make any necessary changes to technical data. If you have technical problems or other questions, please contact our customer service team.

Safety

- The product must not be used by people who are not sensitive to heat or by other very vulnerable people who cannot react to overheating (for example diabetics, people with skin changes caused by medical conditions, people with scars in the treatment area or people who have taken medication or drunk alcohol).
- The product can be used by children aged 8 and over and by people with any form of disability or with a lack of experience and knowledge which could put their safety at risk, provided that they have been given instruction in using the product safely and understand the potential risks. Never let children play with the product. Children must not be allowed to clean or maintain the product unless they are being supervised by an adult.
- Never let children play with the product.
- This product is intended for private domestic use.
- Do not use the product when it is folded or rolled up.
- **Warning:** Do not use the product near baths, showers, washbasins or other containers that hold water or other liquids.
- Consult a doctor if you have a pacemaker or a skin disease or have similar symptoms.
- Never use the product if it is damaged in any way or not working normally.
- The product must only be plugged into to a **5V min. 1A** power source.
- Do not cover the product.
- Do not use the product with other heat sources, such as hot water bottles, heated cushions or similar items.
- Check the product regularly for any signs of wear and tear or damage. If the product is damaged, it must not be used.
- This product is not intended for medical use in hospitals.

Symbols

Keep the product away from sharp objects.

The product must not be used on children under the age of 3.

Do not wash. The product must not be washed.

Do not tumble dry. The product must not be tumble-dried.

Do not iron. The product must not be put through a mangle or ironed.

Do not dry clean. The product must not be dry cleaned.
Do not bleach. Do not use bleach on the product.

Cleaning and maintenance

Note Let the product cool down completely before storing it away or cleaning it.

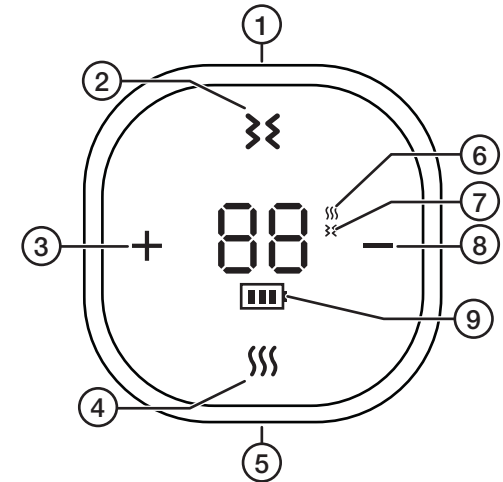
Storage

- Store the product in a dry, dust-free place out of the reach of children if you are not going to use it for a long time.
- Do not put anything on top of the product while it is in storage.

Cleaning

The product will usually only need to be aired, but if it is stained, you can wipe the outside with a damp cloth. Use only a mild detergent. Never use solvents or corrosive chemicals.

Product description

- | | |
|-------------------------|---------------------------------|
| 1. USB-C input port | 6. Heat mode indicator light |
| 2. Massage | 7. Massage mode indicator light |
| 3. + Increase intensity | 8. - Reduce intensity |
| 4. Heat | 9. Battery status indicator |
| 5. Lock/On/Off | |

Charging the product

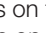
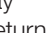
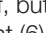
Note Charge the product fully before using it for the first time.

1. Connect a USB-C cable to the input port (1) and to a suitable USB charger.
2. The battery status indicator changes while the battery is charging. The battery status is indicated by the number of static bars. Three static bars means that the battery is fully charged.

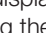
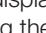
Using the product

Heat

1. Press  to switch the product on. Press  once more to lock the buttons. The display turns off. Press  a third time to unlock the buttons. Press and hold  for 2 seconds to switch the product off.
2. Press  to switch on heat mode. The product starts on the lowest heat. The heat mode indicator light (6) comes on and the display shows **45** (meaning 45°C). After 3 seconds, the display returns to showing the timer.
3. Press  again to switch to medium heat. The display shows **50** (meaning 50°C) for 3 seconds and then returns to showing the timer.
4. Press  again to switch to high heat. The display shows **55** (meaning 55°C) for 3 seconds and then returns to showing the timer.
5. If you press  a fourth time, heat mode switches off, but the product itself does not. The heat mode indicator light (6) goes out.
6. The product has an integrated timer that switches off the heat after 15 minutes.

Massage

1. Press  to switch on the massage function. The product starts in massage mode 1. The massage mode indicator light (7) comes on and **L1** appears on the display. The display returns to showing the timer after 3 seconds.
 - Press **+** or **-** to increase or reduce the intensity.
2. Press  again to switch to massage mode 2. The display shows **L2** for 3 seconds and then returns to showing the timer.
 - Press **+** or **-** to increase or reduce the intensity.
3. Press  again to switch to massage mode 3. The display shows **L3** for 3 seconds and then returns to showing the timer.
 - Press **+** or **-** to increase or reduce the intensity.
4. If you press  a fourth time, massage mode switches off, but the product itself does not. The massage mode indicator light (7) goes out.
5. The product has an integrated timer that switches off the massage function after 15 minutes.

Disposal

This symbol indicates that the product should not be disposed of with general household waste. This applies throughout the entire EU. To prevent any harm to the environment or health hazards caused by incorrect waste disposal, the product must be handed in for recycling so that the materials can be disposed of in a responsible way. To recycle the product, take it to your local recycling facility or contact the retailer. They will ensure that the product is disposed of in a way that does not harm the environment.

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