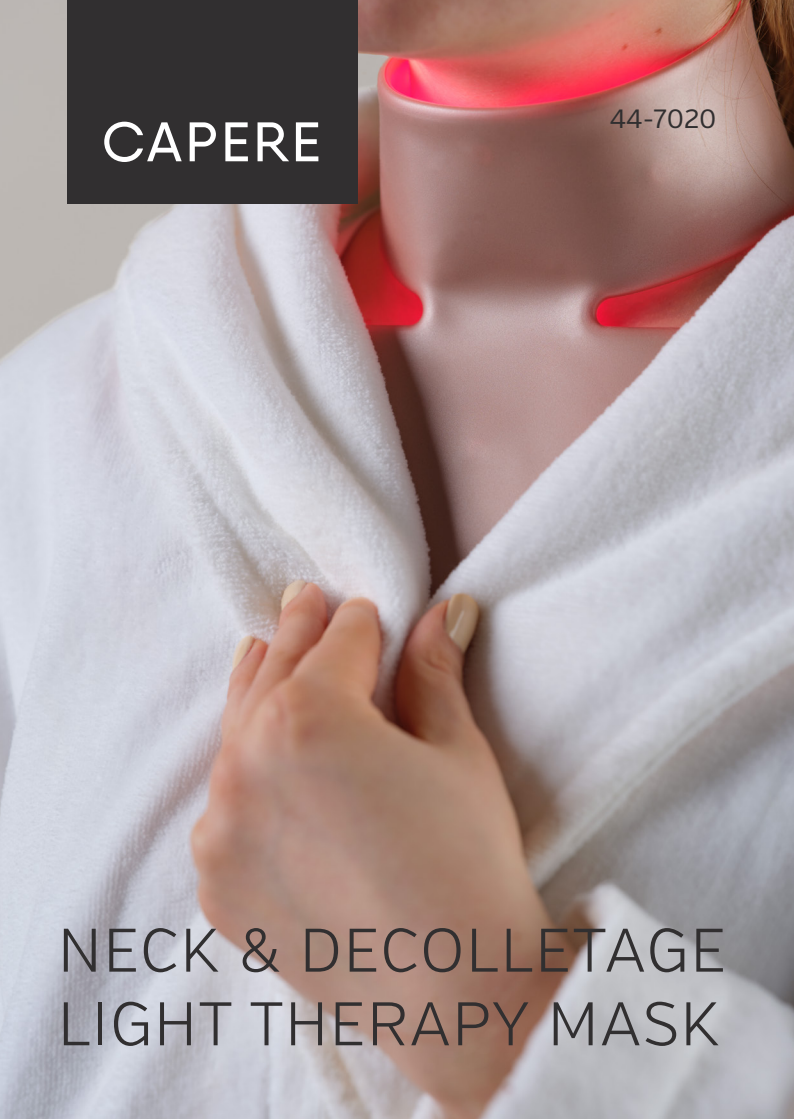
CAPERRE

44-7020

NECK & DECOLLETAGE
LIGHT THERAPY MASK

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Art. no. 44-7020 Model G87

Please read the entire instruction manual before use and save it for future reference. We reserve the right for any errors in text or images and for making any necessary technical changes to this document. If you should have any questions concerning technical problems please contact our Customer Services (see address details on the back).

Safety

- This product is intended for use on the face only.
- The product is not intended for use as a medical treatment.
- The product must not be used by children under the age of 12.
- The product can be used by children aged between 13 and 18 under the supervision of an adult.
- Never use the product in the bathroom or in other damp environments. Never immerse the product in water or other liquid.
- Never use the product if it is malfunctioning or is damaged in any way. Do not use the product if you have dropped it on the floor or if it has fallen into water or other liquid.
- Never stand or sit on the product.
- Never put the product close to a heat source, such as a cooker, woodburning stove, oven or open fire.
- Never cover the product while it is in use.
- Never try to open the housing or repair or modify the product in any way.
- **Risk for eye damage.** Do not look directly at the LEDs.


ENGLISH

- The product is not suitable for use by anyone who has any of the following conditions or who matches any of the following descriptions:
 - Silicone allergy
 - Heart disease
 - Mental illness
 - Abnormal blood pressure
 - Sunburnt skin
 - Keloid scars
 - Oversensitivity to light
 - Discolouration and contact dermatitis
 - Infectious diseases
 - Immunological diseases
 - Diabetes
 - Melanomas on the face
 - Abnormal wound healing
 - Malignant tumours
 - People who are taking medication, cancer patients, especially with skin cancer or precancerous moles
 - People who have had botox injections in the face within the last three months
 - Pregnant women, breastfeeding women and other sensitive people
 - Patients with capillaries caused by long-term use of steroid hormone medications or by liver dysfunction
- If you have permanent or temporary fillers in the area that will be covered, consult a doctor before using the product.
- If you are undergoing hospital treatment, consult a doctor before using the product.
- A certain amount of redness can occur during the treatment. This is normal and should disappear within 1–2 hours. If you feel a burning sensation or other discomfort during the treatment, stop the treatment immediately.

Product description


Control unit

1. USB C connection
2. Battery status indicator
3. Timer indication
4. [NIR] Near infrared (850 nm)
Light mode selector
5. [⌚] Timer, 10, 15, 20 min
6. NIR indicator
7. LED indicator
8. Light mode indicators:
R - Red light (660 nm)
B - Blue light (410 nm)
Y - Yellow light (590 nm)
9. [C] Intensity
10. [⏻] On/Off/Light mode selector

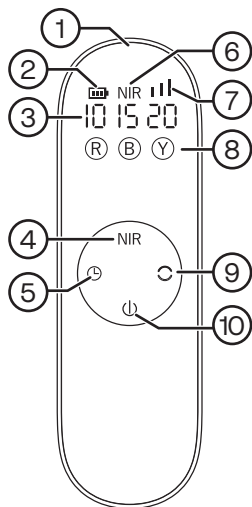Light therapy mask

11. Light therapy mask
12. USB C connection
cable
13. USB A to USB C
connection cable
14. Right strap
15. Left strap

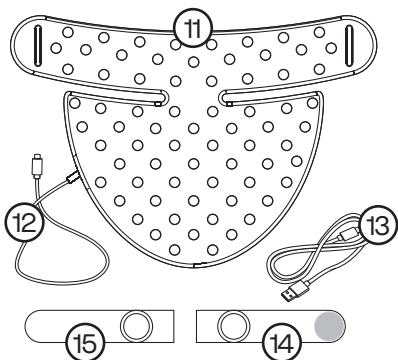Charging

Charge the product fully before using it for the first time.

1. Connect a USB cable to the USB C port (1) on the control box and to a suitable USB charger.
2. The LED indicator (7) wanders while the battery is charging. When the battery is fully charged, only the battery status indicator (2) will be visible
3. The battery status indicator will shine with a red light when the battery needs charging.

Photosensitivity test:

Note You should carry out a photosensitivity test before using the product for the first time:



Switch the product on and put it on the inside of your forearm. Leave it there for 10 minutes on medium intensity (■□□). If any type of redness or rash appears, we advise against using the product, as your skin is probably extra sensitive. Consult a skin specialist before continuing to use the product.

Using the product

Note The product must only be used with the control box provided. The product will be unsafe if you use it in any other way.

Before first use

1. Remove the protective film from both the control unit and the light mask.
2. Attach the straps to the light therapy mask making sure that they are attached to the correct side.
3. Connect the mask's USB-C cable to the USB-C connection on the control unit.

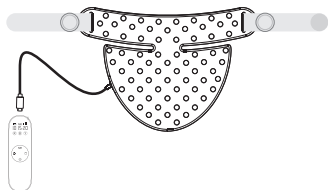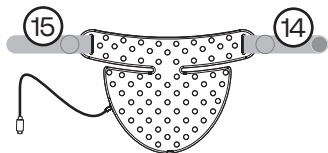Treatment

1. Cleanse your neck and chest. Skincare products or makeup can contain ingredients that have a negative reaction to light.
2. Put the mask on your chest and fasten the straps around your neck with the Velcro fastener. Adjust the mask so that it fits firmly but comfortably.
3. Press and hold [⏻] (10) for 2 seconds to start the mask. The mask starts with red light on the lowest intensity.

ENGLISH

4. Press [] repeatedly to toggle between the following light settings: **Red light (660 nm)** → **Blue light (410 nm)** → **Yellow light (590 nm)** 
5. Press [NIR] to select **Near infrared light (850 nm)**. The NIR indicator lights up. Press [NIR] once more to turn the function off. The NIR indicator goes out.
6. Press [] to select the intensity of the light, the selected intensity is shown with the same number of lit bars in the LED indicator: **Low** () → **Medium** () → **High** () 
7. Press [] to set the duration of the treatment: 10, 15 or 20 min. **Note!** The duration of the blue light treatment should not exceed 10 minutes. To switch off the mask manually, press and hold [] for 2 seconds.
8. You can follow your usual skincare routine after a treatment. **Note** A certain amount of redness can occur during treatment. This is normal and should disappear within 1–2 hours.

Cleaning and maintenance

- Wipe the mask with a soft, slightly damp cloth. Never use solvents, abrasive cleaners or corrosive chemicals. Be careful not to damage the light panel or silicone surface.
- Never immerse the mask in water or rinse it under a tap. It is not waterproof.

Responsible disposal

This symbol indicates that the product should not be disposed of with general household waste.

This applies throughout the entire EU. To prevent any harm to the environment or health hazards caused by incorrect waste disposal, the product must be handed in for recycling so that the materials can be disposed of in a responsible way.

To recycle the product, take it to your local recycling facility or contact the retailer. They will ensure that the product is disposed of in a way that does not harm the environment.

Specifications

Input	5 V DC/1 A
Charging port	USB-C
Charging cable	USB-A – USB-C (included)
Battery	2000 mAh Li-po
Charging time	approx 2 h
Battery life	Low intensity: approx 2 h High intensity: approx 1 h
Wavelength	NIR: 850 nm Red: 660 nm Blue: 410 nm Yellow: 590 nm
Dimensions	Light mask: 347 × 239 × 3.8 mm Control unit: 130 × 40 × 26.5 mm
Weight	Light mask: 193 g Control unit: 91 g

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