CAPERRE

44-9568

LIGHT THERAPY
MASK

Hard light therapy mask

Art. no. 44-9568

Model HML-FM-21V1-7

Please read the entire instruction manual before use and save it for future reference. We reserve the right for any errors in text or images and for making any necessary technical changes to this document. If you should have any questions concerning technical problems please contact our Customer Services (see address details on the back).


Safety

- This product is intended for use on the face only.
- The product is not intended for use as a medical treatment.
- The product must not be used by children under the age of 12.
- The product can be used by children aged between 13 and 18 under the supervision of an adult.
- Never use the product in the bathroom or in other damp environments. Never immerse the product in water or other liquid.
- Never use the product if it is not functioning correctly or is damaged in any way. Do not use the product if you have dropped it on the floor or if it has fallen into water or other liquid.
- Never stand or sit on the product.
- Never put the product close to a heat source, such as a hob, woodburning stove, oven or open fire.
- Never cover the product while it is in use.
- Never try to open the housing or repair or modify the product in any way.
- **Risk for eye damage.** Do not look directly at the LEDs.

ENGLISH

- The product is not suitable for use by anyone with any of the following conditions or who matches any of the following descriptions:
 - Silicone allergy
 - Heart disease
 - Mental illness
 - Abnormal blood pressure
 - Sunburnt skin
 - Keloid scars
 - Oversensitivity to light
 - Discolouration and contact dermatitis
 - Infectious diseases
 - Immunological diseases
 - Diabetes
 - Melanomas on the face
 - Abnormal wound healing
 - Malignant tumours
 - People who are taking medication, cancer patients, especially with skin cancer or precancerous moles
 - People who have had botox injections in the face within the last three months
 - Pregnant women, breastfeeding women and other sensitive people
 - Patients with capillaries caused by long-term use of steroid hormone medications or by liver dysfunction
- If you have permanent or temporary fillers in the area that will be covered, consult a doctor before using the product.
- If you are undergoing hospital treatment, consult a doctor before using the product.
- A certain amount of redness can occur during the treatment. This is normal and should disappear within 1–2 hours. If you feel a burning sensation or other discomfort during the treatment, stop the treatment immediately.

Product description

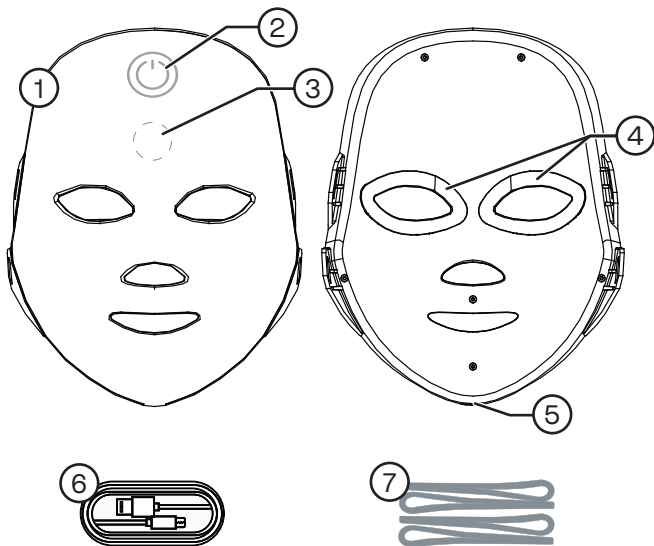1. Light therapy mask
2. Multi-function button
3. LED indicator
4. Eye cups for light protection

5. USB-C input port
6. USB-A to USB-C cable
7. Strap

Charging

Charge the product fully before using it for the first time.

1. Connect a USB cable to the USB-C port (5) on the mask and to a suitable USB charger.
2. The LED indicator (3) lights up red while the product is charging. It changes to blue when the battery is fully charged.

Photosensitivity test

Note You should carry out a photosensitivity test before using the product for the first time:

Switch the product on and put it on the inside of your forearm. Leave it there for 10 minutes. If any type of redness or rash appears, we advise against using the product, as your skin is probably extra sensitive. Consult a skin specialist before continued use.

Using the product

Before using the product for the first time: Attach the strap to the light therapy mask.

1. Cleanse your face. Skincare products or makeup can contain ingredients that have a negative reaction to light.
2. Put the mask on your face and adjust the strap to make sure that it fits firmly but comfortably. **Note:** Make sure that the eye cups fit snugly around your eyes. If they do not, keep your eyes closed during treatment.
3. Press and hold the multi-function button (2) for 3 seconds to start the mask. It will start with red light.

4. Press the multi-function button repeatedly to toggle between 8 different light modes as follows:
Red light → Green light → Blue light → Red + blue light → → Red + green light → Green + blue light → Red + green + blue light → Alternating colours ↻
5. The mask has a timer function that switches it off automatically after 15 minutes. To switch off the mask manually, press and hold the multi-function button for 3 seconds.
6. You can follow your usual skincare routine after a treatment. **Note** A certain amount of redness can occur during treatment. This is normal and should disappear within 1–2 hours.

Cleaning and maintenance

- Clean the mask using a damp cloth. Never use solvents or corrosive chemicals.
- Disinfect the mask regularly to maintain good hygiene. Moisten a cloth with 70% isopropyl alcohol and wipe the inside of the mask. Pay extra attention to the eye cups and other crevices where bacteria can easily accumulate.
- Never immerse the mask in water or rinse it under a tap. It is not waterproof.

Responsible disposal

This symbol indicates that the product should not be disposed of with general household waste.

This applies throughout the entire EU. To prevent any harm to the environment or health hazards caused by incorrect waste disposal, the product must be handed in for recycling so that the materials can be disposed of in a responsible way.

To recycle the product, take it to your local recycling facility or contact the retailer. They will ensure that the product is disposed of in a way that does not harm the environment.

Specifications

Input	5 V DC/2 A
Charging port	USB-C
Charging cable	USB-A – USB-C
Battery	400 mAh
Charging time	approx 1,5 h
Battery life	Single coloured light approx 1,5 h Multi coloured light approx 45 min
Wavelength	Red: 630 nm (± 5 nm) Blue: 460 nm (± 5 nm) Green: 515 nm (± 5 nm) Red/Blue: 630 nm + 460 nm (± 5 nm) Red/Green: 630 nm + 515 nm (± 5 nm) Green/Blue: 515 nm + 460 nm (± 5 nm) Red/Green/Blue: 630 nm + 515 nm + 460 nm (± 5 nm)
Dimensions	232 × 188 × 135 mm
Weight	220 g

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Post Clas Ohlson AB, 793 85 INSJÖN
Osoite Clas Ohlson Oy, Kaivokatu 10 B,